

Video Rhino Laryngoscope

visionflex
Equality healthcare

Ultra-compact laryngoscope with integrated LED and USB.

The Video Rhino Laryngoscope is designed for endoscopy and image capturing by ENT doctors. It is ultra-compact and features an integrated LED light source and USB connectivity.

Available in two different sizes, there is a variant specifically designed for use in paediatrics.

The Video Rhino Laryngoscope was developed in partnership with medical professionals worldwide who have confirmed that the excellent optics, illumination and camera electronics deliver natural and clear images.

Integrate your Video Rhino Laryngoscope with the ProEX Telehealth Hub for the complete telehealth solution.



**Natural and clear
ENT imaging**



Video Rhino Laryngoscope



Ultra-compact laryngoscope with integrated LED and USB.

With precise tactile mechanical probe control, you can navigate the probe tip safely, and intuitively navigate into any position. Smart design and advanced technical construction means less wear and tear and longer product-life.

Razor-sharp images and video capturing with realistic colour reproduction provides the basis for reliable diagnosis and documentation. Capture and record high definition still images and stream

video with the touch of a button on the ProEX, PC, Mac or laptop via USB connection. Images and videos can be archived for later viewing or exported to a network location or USB device. Separate light-source is not required due to the integrated high power LED Advanced Image Technology (AIT).

Available in two different size cameras (3.9mm and 2.9mm diameter) with the 2.9mm specifically designed for use in paediatrics.

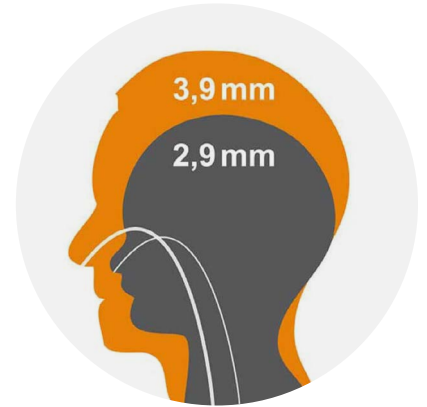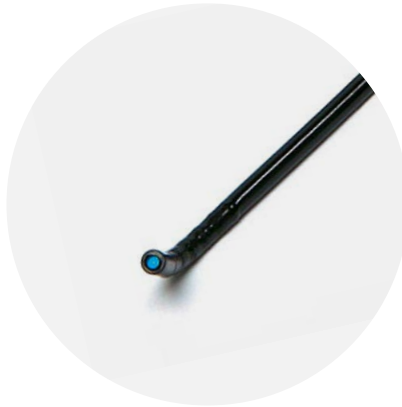

SPECIFICATIONS

Diameter	3.9mm (regular/adult) 2.9mm (small/paediatric)
Field of View (FOV)	90°
Light source	Integrated high power LED lightsource
Tip articulation	130°
Working length	310 mm
Weight	Approx. 330g incl. USB cable
Interface	USB
Cleaning	Completely immersible machine cleaning/disinfection

FEATURES & BENEFITS

- World first ultra-compact scope
- Available in two sizes
- USB connectivity
- Integrated high-power LED
- Advanced Image Technology (AIT)
- Reduced weight and improved mobility
- Compatible with all washing systems and soaking sterilisation, including Tristel Trio50 wipes

For more information or to book a demonstration, please contact Visionflex at sales@visionflex.com

+61 2 8914 4000

visionflex.com



visionflex
Equality healthcare